

DEEP


The PC Deep cushion features a specialized silicone gel insert that molds to the user's ischial tuberosities and coccyx, optimizing pressure redistribution. Its deep-well design with ergonomic contours prevent sacral sitting and enhances stability.

A removable gel overlay sits on top of the base foam delivers immersive comfort.


The moisture-wicking outer cover regulates humidity to maintain skin integrity and minimize pressure injury risk.




Key features / Clinical justification

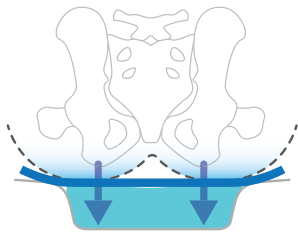

High immersion prevents forward migration of the pelvis

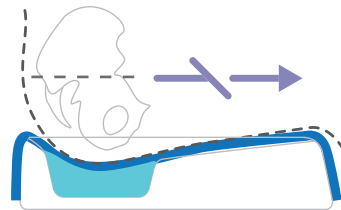

Gel insert in the pelvic well decreases risks of pressure injury on coccyx and sacrum

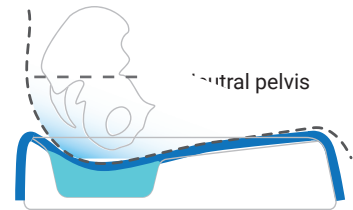

Neutral pelvis and even weight distribution decrease overall risks of pressure injuries

HCPCS E2607, E2608

Custom sizing available upon request

Specification

Sizes	14"W	16"W	18"W	20"W	22"W	24"W
14"D	14x14	16x14	-	-	-	-
16"D	14x16	16x16	18x16	20x16	-	-
18"D	14x18	16x18	18x18	20x18	22x18	24x18
20"D	-	16x20	18x20	20x20	22x20	24x20
22"D	-	-	-	-	22x22	24x22
24"D	-	-	-	-	-	24x24

450lb Max Load